



Product Code

- UTS-0170 Mechanical End-Over-End Shaker
 UTS-0171 Gas Jar with a Rubber Cover for UTS-0170
 UTS-0172 Ground Glass for UTS-0171

| | |
|----------------------------------|------------|
| Models for 220-240V 50 Hz, 1 ph. | UTS-0170-T |
| Models for 110-120V 60 Hz, 1 ph. | UTS-0170-N |
| Models for 220-240V 60 Hz, 1 ph. | UTS-0170-K |

Standards


BS 1377:2

The UTS-0170 Mechanical End-Over-End Shaker is used for the determination of the particle density by the gas jar method and the particle size distribution by sedimentation.

End-Over-End Shaker is used to rotate two gas jars with rubber covers at about 50 rpm.

The 1 L capacity gas jar is made of plexiglass.

UTS-0171 Gas Jar with a rubber cover and a ground glass should be ordered separately.

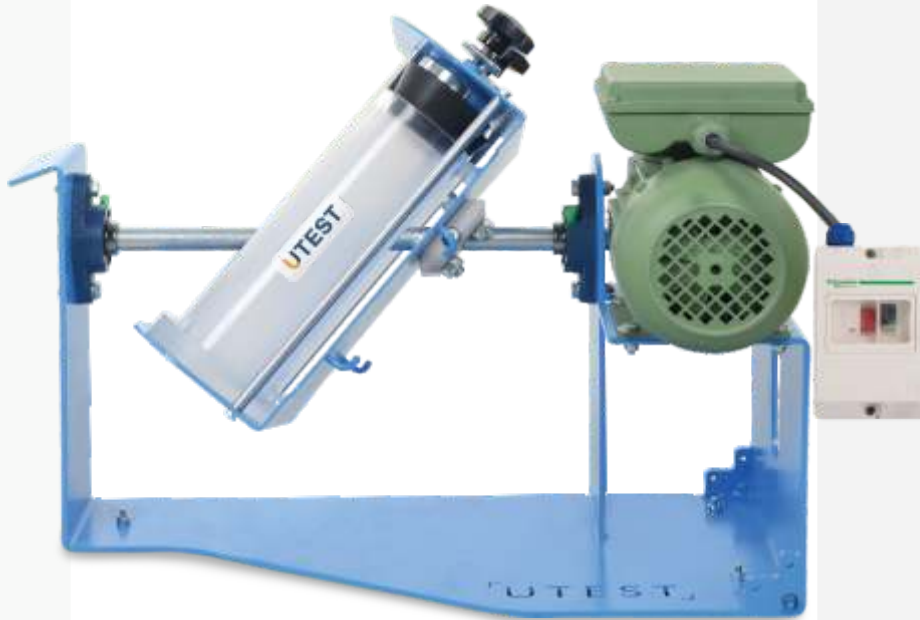

UTS-0170



UTS-0172

UTS-0171

| Product Code | Dimensions | Weight (approx.) | Power |
|--------------|-----------------|------------------|-------|
| UTS-0170 | 500x1030x450 mm | 33 kg | 180 W |
| UTS-0171 | 120x120x280 mm | 0,60 kg | |